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User instructions

Settings before the first use

Please observe the safety regulations

For the first use:

Charge the Relax Smartpillow battery: open the zip on the side of the pillow, locate the cable with micro-USB port and plug in a phone charger, connected to a power supply. Allow 20 minutes for the battery to charge up. (An extension cord is supplied with a universal USB port that can be inserted into any phone jack.) Download and install the Relax Smartpillow

mobile application on your smartphone. The app is available for download on the App Store and Google Play.

Launch the Wifi, your location and the Bluetooth. Make sure you are connected by opening an internet page for example.

Create your account by following the instructions on the smartphone app.

Connect your Relax Smartpillow via Bluetooth on your mobile app. Place a hand on the Relax Smartpillow, it vibrates smoothly to confirm it's connected to your smartphone.

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- Place your Relax Smartpillow on your bed, just as you would a normal pillow.
- When you're ready to go to sleep, open the app and start the recording of the night by clicking on the pillow symbol with green arrow at the right top corner.

- turning the pillow over (this will active Bluetooth signal for 5 minutes).
- Answer the questions on your daily routine
- Set the alarm and the way you prefer to woken up: vibration in the pillow or your phone alarm (this option is only available for Android). You have the option to set a "smart alarm" to wake you a little earlier at the optimum time in your sleep cycle.
- Define a time gap to start the recording later e.g. when you intend to read.
- Click OK to start the monitoring.

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Quick Start Guide 2

- 8. The Relax Smartpillow monitors your sleep and will wake you up at the right time!
- during the night (no radio-waves).
- 10. In the morning, stop the buzzer smart alarm by removing your head from the pillow.
- the Wifi is activated or that you are in 3/4 G. 12. Open the application to end your night by pressing the pillow symbol with the red ung" beschrieben wieder auf. guadrate in the top right corner. The data transfer and synchronisation will begin.

- Approximately after 5 minutes maximum, your night data is available on the application. Let's see how well you slept last night.
- 4. After a week, you will have access to the analysis of your nights.

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Ouick Start Guide 3

- advice on recurring problems).
- Data will be lost.
- your smartphone, additional charges may be incurred for connecting to the Internet or transferring data.

Support

For further product information and frequently asked auestions, please visit our website www. iX21.net. You can also contact us by email on support@ix21.net.

- 5. After 10 nights, the coach will advise you on how to improve the quality of your sleep (a minimum delay is needed to give relevant

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- 5. Check the battery level on the app. There is enough battery power for approximately one week. If the battery level is too low recharge the pillow as described in first use section. Once you have started to monitor a night: do not plug in or recharge the pillow: it would reset the pillow and stop the monitoring.
- **Remark:** Depending on the settings for

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Specifications

- Model: iX-21 Smart Pillow
- iOS 9 and above Wireless compatibility:
- Bluetooth Low Energy (BLE) Wireless range: 5 meters
- Battery life: 5 to 7 days
- Battery: Li-lon battery, rechargeable

Maintenance

Charging instructions Use a USB power supply (e.g. phone charger)

Charging current: 5.0 V

- complying with your country regulations.
- Maximize the battery life by making full battery charges. When the power adapter is inserted into the
- power outlet, do not touch the metal end of the USB plug.

Washing instructions

We recommend spot cleaning with water without soaking the pillow. Your Relax Smartpillow is machine washable at 40°C on a delicate cycle, with low spin drying (< 900 rpm). Do NOT tumble dry.

Disrespecting the instructions of maintenance may damage the product and cause malfunctions.

Software updates: Software updates for the Relax Smartpillow will automatically be provided on your mobile app and installed on the Relax Smartpillow.

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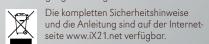
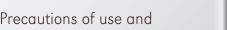
Please note that this pillow device is not a medical device and cannot be used to detect or prevent medical issues. Users with medical illnesses or symptoms should consult a doctor.

Caution: risk of explosion if incorrect power adapter is used. Please follow battery charging instructions (see "Maintenance" section).

Do not use the Relax Smartpillow in an environment with extreme high or low temperatures. Avoid exposing the Relax Smartpillow to direct sunlight. The optimal temperature for the use of this product and accessories is between 16°C and 22°C. The product must not be used under 0°C or over 45°C.

Do not touch the sensor inside the pillow, do not dismantle or try to repair the Relax Smartpillow by yourself, and do not replace the battery.

You should not dispose the pillow in the household waste at the end of the useful life. You are responsible to return used batteries by law to allow professional disposal. In case of questions please contact the local authorities responsible for waste disposal.

Full safety guidelines and instructions can be found at www.iX21.net.

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Legal notice

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that

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The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications or changes to this equipment. Such modifications or changes could void the user's au-

This equipment has been tested and found to

comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular in-

If this equipment does cause harmful interference

thority to operate the equipment.

stallation.

FCC statements

may cause undesired operation.

to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: - Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/ TV technician for help.

EU Manufacturer Declaration

Hereby, Centa-Star declares that this device is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

Any changes or modifications to this equipment not expressly approved by Centa-Star may cause, harmful interference and void the FCC authorization to operate this equipment.